

## Helios 15Q

### 1. Basic product information

Helios 15Q is an ion-exchange resin based on highly cross-linked porous polystyrene-divinylbenzene (PS/DVB) monodisperse microspheres, with a chemically bonded, highly stable hydrophilic layer. Helios 15Q has the characteristics of high flow rate, high resolution, high capacity, high rigidity, good biocompatibility, and physical and chemical stability. Compared with traditional chromatography resin, it significantly improves the production efficiency and sample recovery rate of the downstream purification process, and is widely used in the capture, intermediate purification and fine purification of biomolecules such as antibodies, proteins, peptides and nucleic acids.

### 2. Chromatography resin parameters

|                          |                            |
|--------------------------|----------------------------|
| Resin type               | Strong anion exchange      |
| Matrix                   | Polystyrene-divinylbenzene |
| Median particle size     | 15 µm                      |
| Functional group         | Quaternary amino group     |
| Dynamic binding capacity | ≥45 mg BSA/ml              |
| Recommended flow rate    | 150–900 cm/h               |
| Maximum flow rate        | 1800 cm/h                  |
| Maximum working pressure | 40 bar                     |


### 3. Chemical resistance

|                        |   |
|------------------------|---|
| pH stability*          | 2~12 (working range), 1~14 (CIP)  |
| Chemical stability     | Stable in common aqueous buffers: 1M NaOH+, 1.0 M HCl, 100% ethanol, 100% isopropanol, 100% acetonitrile, 100% methanol |
| Temperature resistance | Don't freeze<br>It can tolerate 121°C high pressure sterilisation   |

\* The physical and chemical properties and functions of the chromatographic resin have no obvious change after being placed in an environment of 40°C and pH 2–12 for 7 days.

### 4. How to use

- Buffer selection: buffer salts whose buffer groups do not act on the resin should be selected. A buffer solution with low salt (less than 5 mS/cm) and low pH (usually 1 pH unit lower than the isoelectric point of the target molecule) should be adopted to facilitate the binding of the target molecule. The stability of samples in the buffer should be considered. Elution buffers are usually made by adding a high concentration of salt (e.g. 1M NaCl) or a low pH buffer to equilibration buffer.
- Flow rate: according to the column bed height, use a flow rate of 150~900 cm/h; the higher the column bed height, the lower the flow rate.
- Sample preparation: to prevent blocking of the column, the sample needs to be filtered by a microporous membrane of 0.45µm before loading. The pH and conductivity of the sample are adjusted to be consistent with the

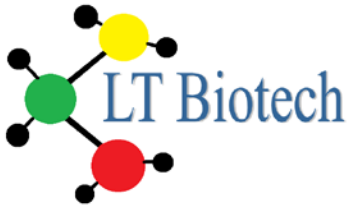

equilibration buffer (the pH and conductivity of the sample can be adjusted by dilution, ultrafiltration and desalination with Persedex G-25).

- Equilibration: wash the column with equilibration buffer until the pH and conductivity of the column outlet buffer are basically the same as the equilibration buffer, which usually needs 3–5 column bed volume (CV).
- Sampling: the loading volume is determined according to the mass of the sample and the binding load of Helios 15Q.
- Rinse: wash the column with equilibration buffer until the UV absorption value is reduced to an appropriate value.
- Elution: a linear gradient or a step-gradient can be used to gradually increase the concentration of salt ions to elute molecules with different binding strengths.
- Regeneration: flush the column with a high concentration of salt (eg: 2M NaCl).
- Re-equilibration: after rinsing with equilibration buffer, the second sample can be loaded and repeated.

## 5. Cleaning and sterilisation

Contaminants (e.g. lipids, endotoxins and proteins) accumulate on the column as the number of uses of the chromatography resin increases. Determine the frequency of CIP according to the degree of contamination of the chromatography resin (if the contamination is considerable, CIP is recommended after each use to ensure repeatability of results and to prolong the working life of the chromatography resin). For different types of impurities and contaminants, the recommended cleaning conditions are as follows:

- Removal of strongly binding proteins: wash with 5 CV of 2M NaCl solution/
- Removal of strongly hydrophobic proteins and precipitated proteins: first wash with 5 CV of 1M NaOH solution, then wash the lye with 5–10 CV of ultra pure or pure water.
- Removal of lipoproteins and lipids: first wash with 5 CV of 70% ethanol or 30% isopropanol, then rinse with 5–10 CV of ultra pure or pure water.

Note: 70% ethanol or 30% isopropanol should be degassed before use; the flow rate should be 30–60 cm/h during CIP. Reverse cleaning should be used when the clogging is severe.

To reduce the microbial load, it is recommended that 0.5~1M NaOH solution is used to treat the chromatographic resin. Treatment time is 15~30 minutes. It can also be autoclaved at 121°C (pH7) for 20 minutes.

## 6. Storage

Keep the unopened chromatography resin in the original container and store at 4~30°C in a well-ventilated, dry and clean place. Do not freeze. Wash the used column with 2–3 CV of 20% ethanol solution and store at 2~8°C.

## 7. Destruction and recycling

Since chromatography resin is difficult to degrade in nature, it is recommended that the waste chromatography resin is incinerated to protect the environment. For chromatography resin that has been in contact with biologically active samples such as viruses and blood, follow the local biosafety requirements before destroying or disposing of it.



LT BIOTECH LTD.  
Mokslininkų g. 6A,  
LT-08412 Vilnius,  
Lithuania  
[www.ltbiotech.lt](http://www.ltbiotech.lt)

## 8. Packing method

Detailed information on resin packaging is available on request. Please contact your local distributor.

## 9. Ordering information

Product name: Helios 15Q

| Product Cat. No | Package |
|-----------------|---------|
| 111-00025       | 25 ml   |
| 111-00100       | 100 ml  |
| 111-00500       | 500 ml  |
| 111-01000       | 1 L     |
| 111-05000       | 5 L     |
| 111-10000       | 10 L    |
| 111-20000       | 20 L    |